

Virtuosity  
is in your hands

**Niti Control,  
the contra-angle with  
adjustable torque**

Located on the dental cart,  
your handpiece is within  
easy reach



**DENTSPLY**  
MAILLEFER

# NiTi Anthogyr Control, maximum efficiency and ergonomic design




## Precision and effectiveness

- 4 torque levels ranging from 0.7 to 4.5 N.cm, for use with all types of Nickel Titanium files
- Automatic declutching if maximum torque is reached

## 2 possible gear ratios

- 1/128 for an electrically operated chair (maximum speed 40'000 t/min)
- 1/64 for a air operated chair (maximum speed 20'000 t/min)

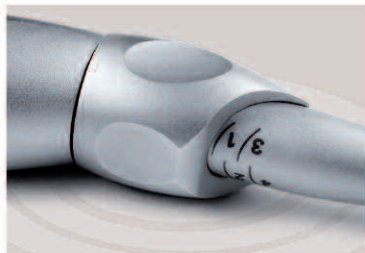

## Ease of use and ergonomic in design

- Automatic push button
- Micro-head: just 9.5 mm for easy access and increased visibility
- A simple setting ring that can be operated with one hand
- Minimum weight: 59 g
- New chrome coating (Corrosion resistance improved)

Distributed by

**DENTSPLY**  
**MAILLEFER**